



DASH FLEX 24 PROGRAMMABLE STAT CENTRIFUGE

The DASH Flex 24 is the solution for high-volume STAT needs. Customize settings and check every detail on the digital display, and save up to 10 custom cycles to produce high-quality STAT samples. Always know when your STAT samples are ready: unique flashing lid lights (patent pending) clearly show cycle status and keep TAT down.

Easy to Monitor

LED lid lighting indicates when the centrifuge is ready (off), running (on), or done (flashing), reducing forgotten samples and keeping your average TAT down. The digital display shows cycle parameters and time remaining.

Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The DASH Flex 24 arrives pre-programmed with the most common STAT settings and can spin 75-100 mm (up to 50 mL) tubes with the same rotor and universal tube holders.

Eliminate Wait Time

The DASH Flex 24 is so compact that it can be used in an array in accessioning. The array will allow spinning samples as soon as they arrive, further shortening your TAT. For a more compact unit, consider the 6-tube DASH Apex 6.

Safe and Durable

Enjoy years of worry-free operation with the DASH's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.

A Great Value

Each centrifuge comes fully equipped and ready to run with a horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.

Tubes	Applications	Centrifugation	Included	Warranty
6, 12 or 24 75-125 mm up to 50 mL	STAT Chemistry, Coag, PPP	Horizontal 3,450 xg Max RCF 4,500 Max RPM	Rotor, Buckets, and Tube Holders	2 years Rotor: lifetime warranty Extended warranty available


Designed, built, and supported in the USA




2-Stage Horizontal Rotor

The DASH Flex 24's patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

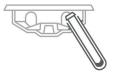

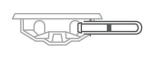

Rotation Angle

DASH FLEX 24

View RPM and remaining spin

time on the digital display



Lid lights indicate cycle status: ready, running, done

Program and name your custom validated cycles Store up to 10 settings

Free 30-Day Evaluation (U.S. Only)

Part Numbers				
Centrifuge	00-184-009-000			
Replacement universal bucket	02-004-0-0012			
1 yr Warranty Extension	03-0-0007-0003			
2 yr Warranty Extension	03-0-0007-0004			
3 yr Warranty Extension	03-0-0007-0005			
* Custom Settings	00-100-101			

* Add this part number to your centrifuge order to have your custom validated cycles pre-programmed at the factory

Specifications				
Sound level	64 dBA			
Dimensions (W x D x H)	14.5 x 17 x 9 in			
	37 x 43 x 23 cm			
Weight	39 lbs (17 kg)			
Power	95 to 253 VAC – 50/60 Hz			

FDA LISTED



RoHS Compliant







druckerdiagnostics.com sales@druckerdiagnostics.com +1-866-265-1486 (U.S. only) +1-814-692-7661

