

Overview



Provide plumbing components, solenoid valve and electrical connections necessary to mix cold tap water with hot drain water so that water $<60^{\circ}\text{C}$ ($<140^{\circ}\text{F}$) is directed to the house drain system. When the washer's drainage pump is activated, the solenoid valve opens releasing the flow of cold tap water to the drain.

Drain Water Cooling Kit provides plumbing components, solenoid valve and electrical connections necessary to reduce drain water temperature to $<140^{\circ}\text{F}$ ($<60^{\circ}\text{C}$). Requires cold water connection. For use in Labconco Glassware Washers manufactured after October 2020. For 230V, 50/60 Hz operation.

Catalog Number: 4678401


Specifications

- **Weight:** 6.0 lbs
- **Weight metric:** 2.7 kg
- **Dimensions:** 3.5" w x 4.0" d x 3.4" h
- **Dimensions metric:** 8.9 x 10.2 x 8.6 cm
- **Electrical:** 208/230V, 50/60 Hz
- **Applicable Products:** FlaskScrubber Glassware Washers, FlaskScrubber Vantage Series Glassware Washers, SteamScrubber Glassware Washers

For Product Assistance



[Amy Ruano](#)
Customer Service Representative

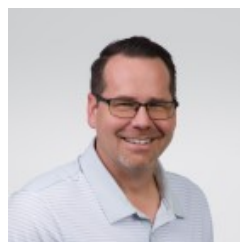

[Dan Simpson](#)
Senior Sales Representative

